

The ultimate protection for Active Pharmaceuticals Ingredients



Protect your valuable sterile API with a proper and RELIABLE primary packaging

Specific Surface treatment

Seamless aluminium cans can be supplied pickled and washed

(ultra-clean washing), or anodised for pharmaceuticals sterile powders as APIs

High Resistance

Seamless aluminium cans are resistant against shocks, perforations, cuts, cold andheat. They also give absolute protection to light and gas permeation.

Quality Standards & Certifications

Metalco is certified ISO 9001:2008, ISO 14001:2004 and EMAS

Quality control

100% on-line inspection of dimensions and surface treatment

Easy sterilisation/ Depyrogenation

- Dry heat sterilisation allows in one single step the sterilisation and depyrogenation of the cans.
- Autoclaving can also be used for sterilisation alone. When tested before dispatch, our cans have been found to be pyrogen free

Guaranteed cleanliness and surface treatment

- Anodization, washing and packing are completely automated
- Final washing and packing is lead in ISO7 clean room
- Surface treatment features are tracked for each container

Batch Unit traceability

Batch and traceability are ensured through ink-jet printing or engraving on the containers

Aluminium recycling and Energy savings

Aluminum is the most recyclable of all materials. Nearly 75 percent of all aluminum produced in the world is still in use today. Aluminum can be recycled directly back into itself over and over again in a true closed loop.

Recycling aluminium uses about 5% of the energy required to create aluminium from bauxite.





ALUminium PACKaging in Pharma

Aluminium canisters are commonly used in pharma sector for:

- β-lactam Antibiotics
- non β-lactam Antibiotics
 Anticancers
- ...and many other sterile APIs.
- Steroid Hormones

Metalco-Alupack provides close regulatory support on aluminium primary packaging.


EASY TO USE CLOSING SYSTEMS

Metalco-Alupack provides accessories in addiction to its containers range. Several closing systems and accessories will bring smooth and reliable use for your confort and safety.



ANODISED PRIMARY PACKAGING

Anodizing is an electrochemical process that converts the Aluminium surface into a durable, corrosion-resistant, anodic oxide finish. The aluminum oxide is not applied to the surface like paint or plating, but is fully integrated with the underlying aluminum substrate, so it cannot chip or peel.

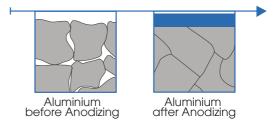



Drug Master File

Metalco-Alupack has registered a Type III (packaging) DMF with the USFDA which is updated and enriched annually.



Accessories and tools for user friendly solutions

Manual closing device

AISI316L stainless steel closing tool allows a high crimping efficiency joined to an entrusting reliability of sealing. Easy to sterilize.





Stopper insert/removal tool

Easy to use device allows to remove and insert the chlorobutyl stopper avoiding any contact with the can's opening part drastically reducing any risk of contamination. Easy to sterilize.

Semi-automatic closing machine

For filling lines Metalco has developed a new concept semi-automatic machine which hermetically crimp and seals the canister.









Metalco products are manufactured according to the strictest american and european specifications and meet the most demanding quality standards

Bureau Veritas ISO 14001



Bureau Veritas ISO 9001



Via Mulino, 44 Montecarlo (LU) Italy

Tel. +39 0583 22090 Tel. +39 0583 229173 Fax +39 0583 22109

info@metalco-alupack.eu www.metalco-alupack.eu